



Key Benefits

- Ultra-reliable connection
- Super simple to use
- Fully monitored
- Long Life Durability



New!

Ultra-Secure EZ CARE Magnetic Socket

EZ CARE's new, revolutionary patient handset magnetic plugs and sockets have been especially developed to set new standards of safety and ease of use. Handset connections to sockets become effortless and super secure. The lifecycle of the devices is also extended.

Handling

The patient handset plug is guided into the socket – the magnetic alignment and clever design ensure automatic, correct positioning. No tilting, no effort required – for maximum ease of use and reliable call triggering.

If the pull is too strong, the plug automatically detaches from the socket without damaging it.

Safety

EZ CARE magnetic sockets provide a reliable, self-adjusting connection that holds firmly and securely through magnetic force. This reduces the risk of short circuits, sparking and mechanical wear.

The plugged-in patient handsets are permanently and fully monitored and meet all necessary safety requirements.

Durability

No worn contacts, no bent plugs – just a reliable, long-lasting solution.

High-quality materials and robust construction guarantee durability, even under demanding conditions.

High reliability – for a future that lasts!

The new plug-in device ...

EZ CARE's innovative plug-in design sets new high standards in healthcare. It delivers maximum peace of mind for staff and patients alike.



The features at a glance:

- Large call button with gender-neutral symbol, red locator light, and reassuring lamp
- Magnetic socket for patient handset devices with and without speech function
- Anti-twist protection thanks to clever design
- Mini-DIN socket for external devices or gray mute button for call mute
- Complete antimicrobial protection
- Fully monitored
- Available in various designs
- Integration into switch ranges possible
- Sustainability through modular systems
- Suitable for wipe disinfection

Reduced risk of infection

The magnetic sockets and patient handsets in the patient room are made exclusively from antimicrobial materials, thereby contributing significantly to reducing the risk of infection.



... with secure connection ...

Reliable, intuitive to use, and fail-safe. For greater safety and comfort for optimal patient care.




Many different designs

EZ CARE offers a wide range of variants with different key combinations and connection options. Depending on requirements, different patient handsets with or without speech functionality are also available for new connectors

Lifecycle friendly, Modular Design

The proven modular system also allows individual components of the new magnetic connectors to be replaced easily and cost-effectively. This means that functional changes can be carried out simply and efficiently during operation.

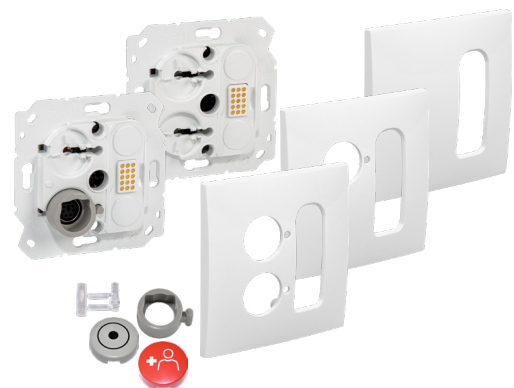

Patient handsets

A wide range of patient handsets with different functions, with or without speech functionality, are available for EZ CARE magnetic sockets. The plug on these patient handheld devices is eject-proof and self-releasing if pulled improperly. For added patient safety reasons, there are no magnets in the plug.



... and proven reliability

EZ CARE stands for a reliable connection between patients and nursing staff. Robust, intuitive, and always ready for use.




Extra protection for patients in intensive care

We also have a dedicated 2 x MOPP interface which a handset can be connected to. The adapter prevents dangerous foreign electrical potentials from being discharged through the patient for extra peace of mind.

Easy repair

Devices damaged through everyday use (e.g., broken button caps or cover plates) can be easily repaired on site by the in-house technician. This saves costs, extends the product life cycle, and thus further increases sustainability.

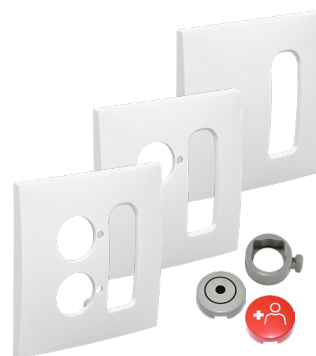

Sustainable and future-proof

The unique backward compatibility of EZ CARE enables the seamless expansion or modernisation of existing installations, during operation and without risk.

With EZ CARE, you protect people and your investments: future-proof, compliant with standards, and economically sustainable.

SCHRACK
SECONET

Schrack Seconet Care Communication GmbH
Eibesbrunnnergasse 18 | 1120 Vienna | Austria

EZ.950.0033E | Version 1.0 | 23. September 2025

© 2025 Schrack Seconet Care Communication GmbH
All rights reserved.

ez-care.com

Company stamp